

WORLD INTELLECTUAL PROPERTY ORGANIZATION  
International Bureau



INTERNATIONAL PATENT CLASSIFICATION <sup>4</sup> : <b>G01N 33/569</b>		<b>A1</b>	(11) International Publication Number: <b>WO 89/ 01162</b> (43) International Publication Date: 9 February 1989 (09.02.89)
(21) International Application Number: <b>PCT/AU88/00274</b> (22) International Filing Date: <b>28 July 1988 (28.07.88)</b> (31) Priority Application Number: <b>PI 3384</b> (32) Priority Date: <b>28 July 1987 (28.07.87)</b> (33) Priority Country: <b>AU</b>		(74) Agent: <b>SPRUSON &amp; FERGUSON; G.P.O. Box 3898, Sydney, NSW 2001 (AU).</b> (81) Designated States: <b>AT (European patent), AU, BE (European patent), CH (European patent), DE (European patent), DK, FI, FR (European patent), GB (European patent), IT (European patent), JP, LU (European patent), NL (European patent), NO, SE (European patent), US.</b>	
(71) Applicant (for all designated States except US): <b>BIO-TECHNOLOGY AUSTRALIA PTY. LTD. [AU/AU]; 28 Barcoo Street, Roseville, NSW 2069 (AU).</b> (72) Inventors; and (75) Inventors/Applicants (for US only) : <b>ATRACHE, Vincent, Habib [AU/AU]; 4 Benaroon Avenue, St. Ives, NSW 2075 (AU). ASH, Megan [AU/AU]; 20/3-7 Edgeworth David Avenue, Hornsby, NSW 2077 (AU). VAN HUYNH, Ca [AU/AU]; 67 Blaxcell Street, Granville, NSW 2142 (AU).</b>		<b>Published</b> <i>With international search report.</i>	
(54) Title: <b>DETECTION METHODS</b>			
(57) Abstract <p>The invention relates to methods for the detection of low numbers of a particular microorganism or microorganisms from a mixed population, the method comprising exposing the sample to a solid support to which antibodies specific for the organism(s) are adsorbed and either growing the bound organism(s) to a detectable level followed by immunoassay or releasing the organisms from the support followed by growth on a non-selective medium and observation of the resultant colonies. The invention also relates to test kits for performing these methods.</p>			
<p align="center"><b>BEST AVAILABLE COPY</b></p>			

***FOR THE PURPOSES OF INFORMATION ONLY***

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	FR	France	ML	Mali
AU	Australia	GA	Gabon	MR	Mauritania
BB	Barbados	GB	United Kingdom	MW	Malawi
BE	Belgium	HU	Hungary	NL	Netherlands
BG	Bulgaria	IT	Italy	NO	Norway
BJ	Benin	JP	Japan	RO	Romania
BR	Brazil	KP	Democratic People's Republic of Korea	SD	Sudan
CF	Central African Republic	KR	Republic of Korea	SE	Sweden
CG	Congo	LI	Liechtenstein	SN	Senegal
CH	Switzerland	LK	Sri Lanka	SU	Soviet Union
CM	Cameroon	LU	Luxembourg	TD	Chad
DE	Germany, Federal Republic of	MC	Monaco	TG	Togo
DK	Denmark	MG	Madagascar	US	United States of America
FI	Finland				

- 1 -

DETECTION METHODSTechnical Field

The current invention relates to methods for detecting low levels of a particular microorganism, or microorganisms from a mixed culture or sample using antibodies and solid immunosorbent supports without the need for a preliminary or further growth step in selective media.

Background Art

Solid-phase immunoassays, based either on enzymes or radioactive isotopes as labels have found wide application in diagnostic microbiology due to their high specificity and sensitivity.

The specificity of an immunoassay is determined by the antibody or antigen which has been immobilized on the solid support. A major advantage of a solid phase assay is that on completion of the immune reaction, unwanted material is easily and rapidly separated from the antigen-antibody immune complex by a simple washing step. A wide variety of solid supports have found application for antibody or antigen immobilization and include polystyrene, polyvinyl chloride, nylon, titanous hydroxide, agarose beads and nitrocellulose.

The successful application of immunoassays to the detection of microorganisms in a sample is possible only if the particular organism of interest is present in sufficient numbers. This critical concentration is determined by the sensitivity of the immunoassay which can vary greatly depending on the affinity and avidity of a particular antibody for its antigen. It is for this reason that many immunoassays require culturing of the sample prior to performing the test.

This usually involves a pre-enrichment step to resuscitate injured microorganisms followed by a selective enrichment to increase the numbers of the microorganism of interest. Selecting culture conditions which favour the growth of a particular microorganism over its competitors has traditionally involved the use of either antibiotics, specific nutritional requirements or manipulation of the physical characteristics of the growth medium, e.g. temperature. These methods can take one to two days or several weeks depending on the organism and the extent of contamination with other microflora.

The use of immunosorbents for selecting microorganisms from a mixed population is known. The sorbent immobilizes a predetermined species of microorganism against which the antibodies are directed. Cells captured in

- 2 -

this way can be incubated in a growth medium and then counted by such traditional techniques as plating and colony counting.

US Patent 4 592 994 describes a method for the determination or identification of microorganisms or unicellular organisms in a sample. The method involves exposing the sample to an adsorbent having "a specific binding power" which may be provided by an antibody raised against the microorganism to be detected. Unbound sample is separated and the adsorbent with bound organisms is exposed to a nutrient medium to initiate metabolism. This nutrient medium undergoes physical or chemical changes as a result of this metabolism and these changes are observed in conjunction with calibration curves to determine presence and amount of the relevant microorganism. The assay involves indirect detection of original numbers of organisms through detection of metabolites in the medium. This may generate problems with regard to the specificity of the assay since different microorganisms may share metabolites.

US Patent 4 563 418 describes a method for the detection of a particular motile organism in a sample, for example, flagellate bacteria such as Salmonella species.

The method involves enriching the sample in an enrichment medium selective for the particular motile organism and filling a motility vessel with a non-selective medium containing a chemotactic attractant which serves to temporarily immobilize the organism of interest and its competitors in the medium for some time after inoculation. Antibodies specific to the flagella of the particular motile microorganism are added through another opening in the motility vessel. The vessel is incubated under sufficient temperature and time conditions to permit the motile organisms to metabolize the chemotactic attractant, thereby reducing its concentration sufficiently to allow movement of the organisms present with the result that the organisms move through the medium and the particular motile organisms being assayed are immobilized by the antibodies. The quantity of antibody used is sufficient to produce a permanent immobilization band.

In Mohit et al ["A Simple Single-step Immunoimmobilization Method for the Detection of Salmonella in the Presence of Large Numbers of Other Bacteria" J. Med. Microbiol. 8 173 (1975)], a method of detecting Salmonella in a mixed population is described. This method employs a selective semi-solid medium, which promotes the migration of Salmonella, followed by immobilization using polyvalent H antisera.

- 3 -

La Roche et al in "Field Evaluation of the Membrane Filter-Disc Immobilization Technique in the Detection of Salmonella in Egg Products" describes a method for detecting Salmonella which involves using a membrane filter to concentrate Salmonella from a primary enrichment broth before selective migration in order to increase recovery.

Stannard reported in the annual report for Leatherhead RA for 1986, investigations into separation of Salmonella from other organisms, with a view to reducing the time required for detection of Salmonella in samples.

This method relied on antibody-coated magnetic particles. However, it was found that the apparent enrichment of Salmonella over other closely related organisms was in fact due to differential affinity of these organisms to the glassware used in the experiment. Attempts to use this differential affinity for glass to enrich Salmonella were unsuccessful due to the lack of specificity of the effect.

Thus it can be seen that the prior art methods for detecting microorganisms from mixed populations provide means for detecting low numbers where enrichment in selective medium is used. Immunoimmobilization has been used but has not been shown to be effective for detecting low numbers in the absence of a such a selection on selective medium.

The current invention provides methods for detecting low numbers of a particular microorganism or microorganisms in a mixed population which overcome the need for pre-selection in selective media by using an immunoimmobilization technique followed by non-selective growth and immunoassay or by cleavage of the antibody-microorganism bond and growth of the microorganisms on non-selective media.

The methods can be especially usefully applied to detecting Salmonella and Listeria spp meeting the need for rapid, sensitive methods for their detection in mixed populations.

#### Description of the Invention

The invention provides methods for the rapid detection of a particular microorganism, or microorganisms, in the presence of competing micro-flora using specific antibodies adsorbed onto a solid support. The immobilized antibody is to an antigen of the microorganism which allows selective capture and immobilization of the desired microorganism without compromising its ability to replicate. Appropriate antibodies can be raised against surface antigens provided by such surface structures as flagella or lipopolysaccharides. Selective concentration of the desired microorganism or microorganisms onto the solid support allows rapid

- 4 -

separation from competing microflora in the sample and is achieved by simply washing the solid support. The immobilized cells can then be transferred to a nutrient broth to allow replication during which time the multiplying microorganisms will continue to be captured and immobilized on the solid support until the antibody sites are saturated.

The time taken for the concentration of the microorganisms to reach a detectable level will depend on the generation time of the particular microorganism or microorganisms and on the sensitivity limits of the immunoassay.

Once a detectable level is reached the solid support is simply separated from the culture broth, then washed and assayed directly.

Alternatively, the antibody-microorganism bond can be broken by an appropriate agent resulting in the release of the microorganism. The released cells can be transferred to a nutrient medium to allow replication and detection of colonies that form.

The selection of the microorganism or microorganisms during the immunoimmobilization step can involve selection of surface antigens that are common to a genus, to a particular species within that genus or if very high selectivity is required, antigens that are specific to a serotype of that species may be chosen.

This technique thus eliminates the need for elaborate and expensive enrichment media in providing rapid sensitive detection of low levels of a particular microorganism or microorganisms from mixed populations.

In a first embodiment, the invention provides a method for detecting low levels of a particular microorganism or microorganisms in the presence of competing microflora in a sample which method comprises: exposing the sample to a solid support to which are adsorbed antibodies specific for the microorganism or microorganisms being detected, said antibodies being capable of selective capture and immobilisation of the microorganism or microorganisms without compromising the ability of the microorganism or microorganisms to replicate; washing the support to remove unbound materials; adding sterile nutrient broth to the support; incubating the support in the nutrient broth at a temperature and for a time relative to the generation time of the microorganism or microorganisms sufficient to allow the microorganism or microorganisms to reach a detectable level; washing the support and then performing an immunoassay on the support using an immunoreagent specific for the microorganism or microorganisms being detected.

- 5 -

Preferably the antibody is raised against a surface antigen of the microorganism. More preferably the surface antigen is a flagellar protein or lipopolysaccharide.

In a preferred form, the antibodies are immobilized onto a solid support comprising a bead, tube or well.

Preferably, the support material is polystyrene, polyvinyl chloride, nylon, titanous hydroxide, agarose beads or nitrocellulose.

Preferred microorganisms according to the invention include Salmonella species and Listeria species.

The exposure of the sample to the antibody-coated support is generally for 1 hour or less. The wash steps are preferably performed using sterile saline buffer, more preferably Tris-saline buffer. Virtually any nutrient broth may be used in the method, however, preferred broths are tryptone soya broth and M broth. The incubation in nutrient medium may be overnight and is usually at 37°C. However, for detection of Salmonella typhimurium, incubation times of approximately 6 hours have been shown to be effective. Where higher numbers of the particular microorganism are present in the original sample incubation times of between 1 and 6 hours can be used. Preferably, the immunoassay is an ELISA.

In a more preferred form the method comprises: exposing a Salmonella test sample for between about 5 and about 30 minutes to a polystyrene support to which are adsorbed anti-Salmonella flagella antibodies, said antibodies being capable of selective capture and immobilisation of Salmonella without compromising the ability of the Salmonella to replicate; washing the support to remove unbound material; adding sterile nutrient broth to the support; incubating the support at 37°C for between about 1 and about 6 hours; washing the support; adding an enzyme-labelled antibody specific to Salmonella; incubating for between about 5 and about 30 minutes; discarding any excess enzyme-labelled antibody; washing the support; adding a chromogenic substrate specific for the enzyme of the enzyme-labelled antibody; and measuring conversion of the chromogenic substrate to a coloured compound.

Preferably the support is a tube, bead or microtitre well.

Preferably the enzyme-labelled antibody is an (anti-Salmonella) antibody - peroxidase conjugate. However, it is recognised that it is possible to use an unlabelled anti-Salmonella antibody in conjunction with a labelled antibody raised against the anti-Salmonella antibody in an indirect version of the assay for bound Salmonella.

- 6 -

Preferably the wash is performed using a Tris-saline solution as wash solution.

Preferably the substrate is an ABTS/H<sub>2</sub>O<sub>2</sub> solution.

Preferably the nutrient broth is tryptone soya broth.

In a second embodiment the invention provides a method for detecting low levels of a particular microorganism or microorganisms in the presence of competing microflora in a sample which method comprises: exposing the sample to a solid support to which are adsorbed antibodies specific for the microorganism or microorganisms, said antibodies being capable of selective capture and immobilisation of the microorganism or microorganisms without compromising the ability of the microorganism or microorganisms to replicate; washing the support to remove unbound material; releasing the bound microorganism or microorganisms with a releasing agent; adding the released material to a nutrient medium; and incubating at a temperature and for a time relative to the generation time of the microorganism or microorganisms being detected to allow the microorganism or microorganisms to reach a detectable level.

Releasing agents suitable for dissociating the antibody-microorganism bond include:

- 1) chaotropic agents such as 4.5 M MgCl<sub>2</sub> pH 7.5 or 2.5M NaI pH 7.5
- 2) polarity reducing agents such as ethylene glycol in solutions of up to 50%
- 3) pH change inducing agents such as glycine/HCl pH 2.5, aqueous NH<sub>3</sub> pH 11 or 0.5% KOH pH 12.5.

Preferably the releasing agent is 0.5% KOH pH12.5. To this solution protein may be added to provide a carrier for the released microorganisms.

In a preferred form the microorganisms are Listeria species or Salmonella species, more preferably Listeria monocytogenes.

It is recognised that prior to the release of the bound microorganisms they may also be subjected to a non-selective growth in nutrient medium to increase their numbers prior to release and detection.

The sample is preferably exposed to the antibody coated support for approximately 5 minutes. The wash solution is preferably a sterile saline solution, more preferably sterile Tris-HCl buffered saline pH7.5. Virtually any nutrient medium may be used, however, Trypticase soya agar is a preferred medium. The plates are preferably incubated at 37° for between 18 and 24 hours. By using nutrient medium in place of selective media incubation can be reduced to as short a time as 12 hours.



- 7 -

In a more preferred form the method comprises: exposing a Listeria test sample for about 5 minutes to a polystyrene support to which are adsorbed antibodies specific to Listeria species, said antibodies being capable of selective capture and immobilisation of Listeria species without compromising the ability of the Listeria species to replicate; washing the support to remove unbound material; releasing the bound Listeria with a releasing agent; adding the released material to a nutrient agar plate; and incubating the plate for about 16-18 hours at 37°C.

Preferably the antibodies are (anti-Listeria flagella) antibodies.

Preferably the releasing agent is a 0.5% KOH solution.

Preferably washing is performed using a Tris/Saline pH 7.5 solution.

Preferably the nutrient agar plate is a trypticase soya agar plate.

The invention also provides a test kit for the detection of low levels of a particular microorganism or microorganisms in a mixed population which kit comprises: a solid support to which are adsorbed antibodies specific for the microorganism or microorganisms being detected; a wash solution; an enzyme-labelled antibody specific for the microorganism or microorganisms being detected; and a solution of a chromogenic substrate for the enzyme of the enzyme-labelled antibody.

In place of the enzyme-labelled antibody the kit may comprise an anti-microorganism antibody together with an enzyme-labelled antibody raised against the anti-microorganism antibody.

The invention also provides a test kit comprising: a solid support to which are adsorbed antibodies specific for the microorganism or microorganisms being detected; a releasing agent; and a wash solution.

#### Best Method of Performing the Invention

The invention is further described with reference to the following examples which are in no way limiting on the scope of the invention.

##### Example 1

##### Selective isolation of Salmonella using antibody coated wells

Salmonella typhimurium (inhouse strain BTA 438) and Citrobacter diversus (BTA 1323) were inoculated into individual M broth solutions and were incubated overnight at 37°C. Each of the cultures grew to  $10^9$  organisms/ml. These were then diluted as follows:

$10^8$  cells/ml C. diversus with  $10^4$  cells/ml S. typhimurium

$10^8$  cells/ml C. diversus with  $10^3$  cells/ml S. typhimurium

- 8 -

- $10^8$  cells/ml C. diversus with  $10^2$  cells/ml S. typhimurium  
 $10^8$  cells/ml C. diversus with  $10^1$  cells/ml S. typhimurium  
 $10^8$  cells/ml C. diversus with  $10^0$  cells/ml S. typhimurium  
 $10^8$  cells/ml C. diversus only (control)

Polystyrene wells coated with highly purified antibodies to Salmonella flagella were then added to the M broths containing the mixed cultures. Wells coated with non-immune sheep immunoglobulin were placed into a duplicate set of mixed cultures as a control. The broths containing the wells were shaken for 1 hour at 37°C.

The wells were then carefully removed and washed 3 times in Tris-saline, then placed in fresh M broth solutions. These M broths were incubated overnight at 37°C to allow the immobilised organisms to multiply. The wells were then removed from the M broths and an ELISA was performed on them.

#### EIA Method

The wells were washed three times with Tris-saline-Tween (TST), then anti-Salmonella IgG - Horseradish Peroxidase labelled conjugate was added. After 30 minutes incubation at 37°C, the wells were washed three times with TST, then substrate (2,2-azinobis(3-ethylbenzthiazoline sulfonic acid)) in citrate-phosphate buffer (0.1M, pH4) containing hydrogen peroxide (0.005%) was added and colour allowed to develop.

ORGANISMS/ML			OPTICAL DENSITY <sup>1</sup>	
C. diversus		S. typhimurium	<u>Salmonella</u> Ab Coated Wells	Control Wells <sup>2</sup>
$10^8$	+	$10^4$	2.0	0.026
$10^8$	+	$10^3$	1.378	0.035
$10^8$	+	$10^2$	0.306	0.063
$10^8$	+	$10^1$	0.314	0.042
$10^8$	+	$10^0$	0.066	0.050
$10^8$			0.078	0.037

- Optical densities read using dual wavelength ELISA reader at 414nm and 490nm.
- Non-immune sheep immunoglobulin coated wells.

#### Example 2

The method was compared with the Standard Culture Enrichment protocol [AOAC Official Methods of Analysis 963-971 (1984)] and showed increased speed of performance, and enhanced selectivity and sensitivity as described below.

- 9 -

The protocol was as described in Example 1 except that 1ml of the M broth containing the mixed cultures was placed into 10ml of tetrathionate broth and incubated at 37°C for 6 hours.

Following this, an ELISA was performed on the broths. To the remaining 9ml of broth, wells coated with highly purified antibodies to Salmonella flagella were added and the procedure of Example 1 followed except that the incubation of the M broth was for 6 hours only and then an ELISA was performed.

ORGANISMS/ML			OPTICAL DENSITY <sup>1</sup>	
C.DIVERSUS		S.TYPHIMURIUM	SALMONELLA Ab COATED WELLS	TETRATHIONATE BROTHS
10 <sup>8</sup>	+	10 <sup>4</sup>	2.0	0.113
10 <sup>8</sup>	+	10 <sup>3</sup>	2.0	0.089
10 <sup>8</sup>	+	10 <sup>2</sup>	0.244	0.076
10 <sup>8</sup>	+	10 <sup>1</sup>	0.200	0.061
10 <sup>8</sup>		control	0.030	0.065

<sup>1</sup>Optical densities read using dual wavelength ELISA reader at 414nm and 490nm.

After 6 hours the ELISA detected 10<sup>3</sup> Salmonella in the presence of 10<sup>8</sup> Citrobacter from the coated wells whilst after 6 hours selective enrichment in tetrathionate no Salmonella were detected by the ELISA.

### Example 3

#### Selective isolation of Listeria monocytogenes from mixed cultures

Listeria monocytogenes (BTA No. 1767), Staphylococcus aureus (BTA No. 1414), and Streptococcus faecalis (ATCC No. 19433) were individually cultured in Tryptose-Soya broth overnight at 28°C. Each of the cultures grew to a minimum of 10<sup>9</sup> organisms/ml and were then diluted and mixed as follows:

Mixture A: 10<sup>9</sup> cells/ml Staphylococcus aureus  
+ 10<sup>6</sup> cells/ml L. monocytogenes.

Mixture B: 10<sup>9</sup> cells/ml Streptococcus faecalis  
+ 10<sup>6</sup> cells/ml L. monocytogenes.

Controls: 10<sup>9</sup> cells/ml S. faecalis  
10<sup>9</sup> cells/ml S. aureus  
10<sup>6</sup> cells/ml L. monocytogenes

A 200μl sample from each of the above cultures was incubated in a

- 10 -

polystyrene well, precoated with highly purified antibodies to Listeria flagella. After 5 minutes incubation at room temperature the wells were emptied and washed three times with sterile Tris-HCl buffered saline pH 7.5. Any captured organisms were then released into solution by the addition of 0.5% KOH solution pH 12.5 (100µl) and immediately transferred to Tryptose Soya Agar plates and spread using a glass spreader. For comparative purposes each of the cultures above was also examined directly without immunoenrichment by streaking a 10µl aliquot onto Tryptose Soya Agar. All plates were incubated for 18-24 hours at 37°C and then examined under oblique light for typical colonies of Listeria monocytogenes (sparkling blueish-grey, translucent, 0.5-1.5mm in diameter, watery consistency).

### RESULTS

Sample Tested	Colonies Counted			
	Immuno-enrichment		Direct Plating	
	<u>Listeria</u>	Other	<u>Listeria</u>	Other
Mixture A	253	222	0	>1000
Mixture B	125	264	0	>1000
<u>S. aureus</u> 10 <sup>9</sup> cells/ml	-	246	-	>1000
<u>S. faecalis</u> 10 <sup>9</sup> cells/ml	-	500	-	>1000
<u>L. monocytogenes</u> 10 <sup>6</sup> cells/ml	127	-	>1000	-

In a single step, approximately a one thousand fold enrichment of L. monocytogenes over the competing microorganism was achieved. Isolation was greatly simplified and detection was possible when direct plating of such a mixture failed to isolate a Listeria colony.

#### Example 4

#### Sensitivity of Immunoenrichment for the isolation of Listeria monocytogenes

Listeria monocytogenes was cultured in Tryptose Soya Broth for 18 hours at 28°C. The culture was then diluted in sterile saline to give concentrations of 10<sup>7</sup>, 10<sup>6</sup>, 10<sup>5</sup> and 10<sup>4</sup> organisms/ml. Samples (200µl) of each dilution, tested in triplicate, were incubated in individual anti-Listeria antibody coated wells for 5 minutes then the wells emptied and washed with sterile buffered saline pH 7.5. Any immobilized microorganisms were released into solution by the addition of 0.5% KOH solution pH 12.5 and then immediately transferred to Tryptose Soya Agar

- 11 -

plates and spread using a glass spreader. Plates were incubated at 37°C for 18-24 hours and the number of colonies counted.

### RESULTS

<u>Number of organisms added/well</u>	<u>No. of Colonies counted</u>
$2 \times 10^6$	1000
$2 \times 10^5$	300
$2 \times 10^4$	50
$2 \times 10^3$	1

Sensitivity of the method allowed detection of Listeria monocytogenes when as few as  $2 \times 10^3$  organisms were incubated in the antibody coated well.

### Example 5

#### General Method:

A Listeria test sample was adjusted to pH 7.5-8.0 and 200µl were transferred into a USA Dynatech Removawell™ which has been coated with anti-Listeria antibodies in 10mm phosphate buffer pH 8.1, 10µg/ml O/N at 20-25°C.

The sample was incubated at room temperature for 5 minutes.

The well was then washed gently with sterile Tris-Saline solution at pH 7.5.

100µl of KOH (pH 11.1) was added.

A 100µl aliquot of the resulting KOH mix was transferred to a Tryptone soya agar or modified McBride agar plate and spread using a glass spreader. The plate was incubated for between 18 and 24 hours at 37°C and the plate examined for typical Listeria colonies.

#### Optimization of reaction time for Immunoenrichment

Culture: Listeria monocytogenes (BTA 1767) was grown in OXOID's TSB + 0.6% yeast extract at 28°C overnight.

Releasing reagent: 0.5% KOH (pH 10.5).

Plate: Modified McBride Agar plate.

- 12 -

Microorganisms	No. of <i>Listeria</i> colonies isolated on MMA plate				
	5 min.	10 min.	15 min.	20 min	30 min.
* $10^5$ cells/ml of <i>L. m.</i>	58	54	37	7	15
* $10^5$ cells/ml of <i>L. m.</i> + $10^9$ cell/ml					
<i>Streptococcus faecalis</i>	100	110	84	-	10

Immuno-enrichment for *Listeria* using different releasing reagents.

<i>Listeria</i> monocytogenes		No. of colonies on TSA plate
dilution	Releasing reagent	
( $2 \times 10^6$ cell/ml)	0.01% KOH + 0.85% NaCl	145
"	" "	149
<hr/>		
( $2 \times 10^6$ cell/ml)	0.01% NaOH + 0.85% NaCl	191
"	" "	260
<hr/>		
( $10^6$ cell/ml)	0.01% KOH + 0.85% NaCl	122
<hr/>		
"	0.01% NaOH + 0.85% NaCl	145
"	" "	118
<hr/>		
"	10mM Tris (pH 11.1) + 0.85% NaCl	112
"	" "	161
<hr/>		
<i>Listeria innocua</i> (CT-94)		No. of colonies on TSA plates
	Releasing reagent	
"	0.01% KOH + 0.85% NaCl	58
"	" "	65
<hr/>		
( $10^6$ cell/ml)	0.1M Na <sub>2</sub> CO <sub>3</sub> (pH 11.1)	163
"	" "	126

- 13 -

Specificity of Immunoenrichment(1) Culture: All cultures were grown in OXOID's TSB at 28°C overnight.Releasing reagent: 0.01% KOH + 0.85% NaCl (pH 11.1)Microtitre wells: Listeria antibody coated wells.Staphylococcus enterotoxin antibody coated well.

Blank well (normal gamma irradiated polystyrene well).

Mixture of micro-organisms	Type of well	No. of Listeria on TSA plate
* $10^8$ cell/ml <u>L. monocytogenes</u>	<u>Listeria</u> antibody	2,000-3,000
+ $10^8$ cell/ml of <u>Rhodococcus</u>	coated well	
<u>equi</u> + $10^9$ cell/ml <u>Staphylococcus</u>		2,000-3,000
* The same mixture as above	Staphylococcus extero-toxin antibody coated well.	0
		0
* The same mixture as above	Blank well	0
		1

Example 6Salmonella Rapid Detection Kit

This kit enables the user to rapidly test for the presence of Salmonella in food, environmental or clinical samples. Highly specific antibodies to Salmonella flagella coated onto the surface of polystyrene tubes are used to capture Salmonella organisms present in a test sample. This is achieved by incubating the test sample in the antibody coated tube for a short time e.g: 5-30 minutes. The tube is then washed thoroughly to remove unbound material. Sterile nutrient broth e.g: Tryptone soya broth is then added and the tube incubated at 37°C for 1-6h. During this time the immobilised Salmonella replicate and continue to be captured by available antibodies on the tube surface. This allows for a sufficient concentration of organisms to be reached for subsequent detection by immunoassay.

The presence of the captured organisms can now be detected by discarding the culture broth and then adding an enzyme-labelled antibody specific to Salmonella. After a short incubation (5-30 min) excess

- 14 -

enzyme-antibody reagent is discarded and the tube washed.

The enzyme-antibody conjugate which has specifically bound to the immobilised organisms is detected by addition of a substrate for the enzyme.

#### Materials Provided

- anti-Salmonella antibody coated polystyrene tubes
- anti-Salmonella-antibody-Peroxidase conjugate
- wash solution Tris Saline Tween
- Substrate solution - ABTS/H<sub>2</sub>O<sub>2</sub>

By use of Listeria reagents in the place of the Salmonella reagents this kit can be adapted to the detection of Listeria.

#### Listeria Detection Kit

The kit allows the selective isolation of Listeria from culture broths derived from either food or environmental samples and eliminates the need to develop/use selective agar plates for the isolation of the organism.

Highly specific antibodies to Listeria are immobilised on the internal surface of polystyrene tubes provided. An aliquot of the test solution is incubated in the antibody coated tube for a specified time (5 minutes). The tube is then emptied and washed to remove unbound material. Any captured Listeria is then released by the addition of a releasing agent e.g. 0.5% KOH solution. The Listeria released into solution are spread onto a nutrient agar e.g. Trypticase soya agar and the plate incubated for 16-18 hours at 37°C to allow the organisms to replicate.

By this procedure the desired organism is readily isolated from an initial culture broth, (sample) contaminated with a variety of other organisms and thus obviates the need for selective agars. In the specific case of Listeria the described method shows superior selectivity over the selective agars currently recommended by the Microbiological Standard Methods.

#### Materials Provided

- (1) Anti-Listeria flagella antibody coated polystyrene tubes
- (2) Releasing agent (0.5% KOH solution)
- (3) Wash Buffer Tris/Saline pH 7.5

By replacing the Listeria reagents with Salmonella reagents this kit can be adapted to the detection of Salmonella.

These examples illustrate the usefulness of these methods for the rapid detection of a specific organism in a mixed culture.



- 15 -

The methods are more rapid than traditional culture enrichment protocols, with a positive result being obtainable as rapidly as 6 hours, compared to a minimum of 48 hours by standard culture methods.

Industrial Application

The current invention provides an alternative to use of specialized selection media for the detection of low numbers of a particular microorganism or microorganisms in a mixed population.

- 16 -

CLAIMSDETECTION METHODS

1. A method for the detection of low levels of a particular microorganism or microorganisms in the presence of competing microflora in a sample which method comprises: exposing the sample to a solid support to which are adsorbed antibodies specific for the microorganism or microorganisms being detected, said antibodies being capable of selective capture and immobilisation of the microorganism or microorganisms without compromising the ability of the microorganism or microorganisms to replicate; washing the support to remove unbound materials; adding sterile nutrient broth to the support; incubating the support in the nutrient broth at a temperature and for a time, relative to the generation time of the microorganism or microorganisms, sufficient to allow the microorganism or microorganisms to reach a detectable level; washing the support; and then performing an immunoassay on the support using an immuno-reagent specific for the microorganism or microorganisms being detected.
2. The method according to claim 1, wherein the antibodies are antibodies raised against a surface antigen of the microorganism or microorganisms.
3. The method according to claim 2, wherein the surface antigen is a flagellar protein or lipopolysaccharide.
4. The method according to any one of claims 1 to 3, wherein the solid support is a bead, tube or well.
5. The method according to claim 4, wherein the support is polystyrene, polyvinyl chloride, nylon, titanous hydroxide, agarose beads or nitrocellulose.
6. The method according to any one of claims 1 to 5, wherein the microorganism is a Salmonella species or a Listeria species.
7. The method according to any one of claims 1 to 6, wherein the sample is exposed to the support for one hour or less.
8. The method according to any one of claims 1 to 7, wherein the support is washed with a sterile saline buffer.
9. The method according to claim 8, wherein the sterile saline buffer is Tris-Saline buffer.
10. The method according to any one of claims 1 to 9, wherein the nutrient broth is tryptone soya broth or M broth.
11. The method according to any one of claims 1 to 10, wherein the incubation in nutrient broth is overnight at 37°C.

- 17 -

12. The method according to any one of claims 1 to 11, wherein the incubation in nutrient broth is at 37°C for 6 hours.

13. The method according to any one of claims 1 to 12, wherein the immunoassay is an ELISA.

14. The method according to claim 13, wherein the immunoreagent is an enzyme-labelled anti-microorganism antibody.

15. The method according to claim 13, wherein the immunoreagent comprises an anti-microorganism antibody and an enzyme-labelled antibody raised against the anti-microorganism antibody.

16. The method according to any one of claims 1 to 14 which method comprises: exposing a Salmonella test sample for between about 5 and about 30 minutes to a polystyrene support, to which are adsorbed anti-Salmonella flagella antibodies, said antibodies being capable of selective capture and immobilisation of Salmonella without compromising the ability of the Salmonella to replicate; washing the support to remove unbound material with a first wash solution; adding sterile nutrient broth to the support; incubating the support at 37°C for between about 1 and about 6 hours; washing the support with a second wash solution; adding an enzyme-labelled antibody specific to Salmonella; incubating for between about 5 and about 30 minutes; discarding any excess enzyme-labelled antibody; washing the support with the second wash solution; adding a chromogenic substrate specific for the enzyme of the enzyme-labelled antibody; and measuring conversion of the chromogenic substrate to a coloured compound.

17. The method according to claim 16, wherein the support is a tube, bead or microtitre well.

18. The method according to claim 16 or claim 17, wherein the enzyme-labelled antibody is an (anti-Salmonella) antibody-peroxidase conjugate.

19. The method according to any one of claims 16 to 18, wherein the first wash is performed using a Tris-Saline solution as wash solution and the second wash is performed using a Tris-Saline-Tween solution as wash solution.

20. The method according to any one of claims 16 to 19, wherein the substrate is an ABTS/H<sub>2</sub>O<sub>2</sub> solution.

21. The method according to any one of claims 16 to 20, wherein the nutrient broth is tryptone soya broth.

22. A method for the detection of low levels of a particular microorganisms or microorganisms in the presence of competing microflora in

- 18 -

a sample which method comprises: exposing the sample to a solid support to which are adsorbed antibodies specific for the microorganism or microorganisms, said antibodies being capable of selective capture and immobilisation of the microorganism or microorganisms without compromising the ability of the microorganism or microorganisms to replicate; washing the support to remove unbound material; releasing the bound microorganism or microorganisms with a releasing agent; adding the released material to a nutrient medium; and incubating at a temperature and for a time relative to the generation time of the microorganism or microorganisms being detected to allow the microorganism or microorganisms to reach a detectable level.

23. The method according to claim 22, wherein the releasing agent is a chaotropic agent, polarity reducing agent or pH change inducing agent.

24. The method according to claim 23, wherein the releasing agent is a chaotropic agent selected from 4.5 M  $\text{MgCl}_2$  pH 7.5 or 2.5M NaI pH 7.5 or a polarity reducing agent selected from ethylene glycol in solutions of up to 50% or a pH change inducing agent selected from glycine/HCl pH 2.5, aqueous  $\text{NH}_3$  pH 11 or 0.5% KOH pH 12.5.

25. The method according to any one of claims 22 to 24, wherein the microorganism is a Listeria species or Salmonella species.

26. The method according to claim 25, wherein the Listeria species is Listeria monocytogenes.

27. The method according to any one of claims 22 to 26, wherein the sample is exposed to the support for approximately 5 minutes.

28. The method according to any one of claims 22 to 27, wherein the support is washed with sterile buffered saline.

29. The method according to claim 28, wherein the sterile buffered saline is Tris-HCl buffered saline pH 7.5.

30. The method according to any one of claims 22 to 29, wherein the nutrient medium is Trypticase soya agar.

31. The method according to any one of claims 22 to 30, wherein the nutrient medium is incubated at 37°C for between about 12 and about 24 hours.

32. The method according to any one of claims 22 to 31, wherein the nutrient medium is incubated at 37°C for between about 18 and about 24 hours.

33. The method according to any one of claims 22 to 32, wherein the method comprises: exposing a Listeria test sample for about 5 minutes to a polystyrene support to which are adsorbed antibodies specific to Listeria

- 19 -

species, said antibodies being capable of selective capture and immobilisation of Listeria species without compromising the ability of the Listeria species to replicate; washing the support to remove unbound material; releasing the bound Listeria with a releasing agent; adding the released material to a nutrient agar plate; and incubating the plate for about 16-18 hours at 37°C.

34. The method according to claim 33, wherein the nutrient agar plate is a trypticase soya agar plate.

35. The method according to claim 33 or claim 34, wherein the antibodies are (anti-Listeria flagella) antibodies.

36. The method according to any one of claims 33 to 35, wherein the releasing agent is a 0.5% KOH pH 12.5 solution.

37. The method according to any one of claims 33 to 36, wherein the washing is performed using a Tris/Saline pH 7.5 solution.

38. A test kit for the detection of low levels of a particular microorganism or microorganisms in a mixed population which kit comprises: a solid support to which are adsorbed antibodies specific for the microorganism or microorganisms being detected, said antibodies being capable of selective capture and immobilisation of the microorganism or microorganisms without compromising the ability of the microorganism or microorganisms to replicate; an enzyme-labelled antibody specific for the microorganism or microorganisms being detected; and a solution of a chromogenic substrate for the enzyme of the enzyme-labelled antibody.

39. The test kit according to claim 38, wherein the support comprises anti-Salmonella antibody coated polystyrene tubes.

40. The test kit according to claim 38 or claim 39, wherein the enzyme-labelled antibody comprises anti-Salmonella-antibody-Peroxidase conjugate.

41. The test kit according to claims 38 to 40, wherein the substrate solution comprises ABTS/H<sub>2</sub>O<sub>2</sub>.

42. A test kit for the detection of low levels of a particular microorganism or microorganisms in a mixed population which kit comprises: a solid support to which are adsorbed antibodies specific for the microorganism or microorganisms being detected, said antibodies being capable of selective capture and immobilisation of the microorganism or microorganisms without compromising the ability of the microorganism or microorganisms to replicate; an antibody specific for the microorganism or microorganisms; an enzyme-labelled antibody specific for the antibody

- 20 -

specific for the microorganism or microorganisms; and a solution of a chromogenic substrate for the enzyme of the enzyme-labelled antibody.

43. A test kit comprising: a solid support to which are adsorbed antibodies specific for the microorganism or microorganisms being detected, said antibodies being capable of selective capture and immobilisation of the microorganism or microorganisms without compromising the ability of the microorganism or microorganisms to replicate; a releasing agent and a wash solution.

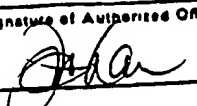
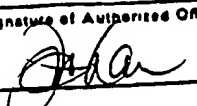
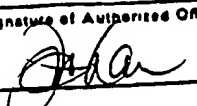
44. The test kit according to claim 43, wherein the support comprises anti-Listeria flagella antibody coated polystyrene tubes.

45. The test kit according to claim 43 or claim 44, wherein the releasing agent is 0.5% KOH solution.

46. The test kit according to claim 45, wherein the wash solution comprises Tris/Saline pH 7.5.

# INTERNATIONAL SEARCH REPORT

International Application No PCT/AU 88/00274

<b>I. CLASSIFICATION OF SUBJECT MATTER</b> : 1 search classification symbols apply, indicate all According to International Patent Classification (IPC) or to both National Classification and IPC Int. Cl. <sup>4</sup> G01N 33/569																							
<b>II. FIELDS SEARCHED</b> Minimum Documentation Searched * Classification System Classification Symbols CHEMICAL ABSTRACTS WPI, WPIL, USPA : Keywords : Microorganisms or (Derwent Databases) Salmonella or Listeria and Immunoassay Documentation Searched other than Minimum Documentation to the extent that such Documents are included in the Fields Searched * AU : IPC G01N 33/569, 33/54																							
<b>III. DOCUMENTS CONSIDERED TO BE RELEVANT</b> <table border="1"> <thead> <tr> <th>Category *</th> <th>Citation of Document, ** with indication, where appropriate, of the relevant passages **</th> <th>Relevant to Claim No. **</th> </tr> </thead> <tbody> <tr> <td>X</td> <td>AU,A, 53187/86 (UNISEARCH LIMITED)(IBRAHIM, George et al) 31 July 1986 (31.07.86)</td> <td>(1-6,16-17, 38)</td> </tr> <tr> <td>P,A</td> <td>US,A, 4752562 (SHEIMAN et al) 21 June 1988 (21.06.88)</td> <td>(1,4-5,22,38, 42)</td> </tr> <tr> <td>A</td> <td>US,A, 4434236 (FREYTAG) 28 February 1984 (28.02.84)</td> <td>(1,5)</td> </tr> <tr> <td>A</td> <td>WO,A, 84/02193 (E.I. DUPONT DE NEMOURS AND COMPANY) (EBERSOLE et al) 7 June 1984 (07.06.84)</td> <td>(1-2,4-5)</td> </tr> <tr> <td>A</td> <td>AU,A, 66596/86 (ORTHODIAGNOSTIC SYSTEMS INC)(FINKENAU) 18 June 1987 (18.06.87)</td> <td>(1-2,4-5,13, 38)</td> </tr> <tr> <td>A</td> <td>EP,A, 0201211 (WHITTAKER CORPORATION)(KHEILL et al) 12 November 1986 (12.11.86)</td> <td>(1,4-6,22)</td> </tr> </tbody> </table>			Category *	Citation of Document, ** with indication, where appropriate, of the relevant passages **	Relevant to Claim No. **	X	AU,A, 53187/86 (UNISEARCH LIMITED)(IBRAHIM, George et al) 31 July 1986 (31.07.86)	(1-6,16-17, 38)	P,A	US,A, 4752562 (SHEIMAN et al) 21 June 1988 (21.06.88)	(1,4-5,22,38, 42)	A	US,A, 4434236 (FREYTAG) 28 February 1984 (28.02.84)	(1,5)	A	WO,A, 84/02193 (E.I. DUPONT DE NEMOURS AND COMPANY) (EBERSOLE et al) 7 June 1984 (07.06.84)	(1-2,4-5)	A	AU,A, 66596/86 (ORTHODIAGNOSTIC SYSTEMS INC)(FINKENAU) 18 June 1987 (18.06.87)	(1-2,4-5,13, 38)	A	EP,A, 0201211 (WHITTAKER CORPORATION)(KHEILL et al) 12 November 1986 (12.11.86)	(1,4-6,22)
Category *	Citation of Document, ** with indication, where appropriate, of the relevant passages **	Relevant to Claim No. **																					
X	AU,A, 53187/86 (UNISEARCH LIMITED)(IBRAHIM, George et al) 31 July 1986 (31.07.86)	(1-6,16-17, 38)																					
P,A	US,A, 4752562 (SHEIMAN et al) 21 June 1988 (21.06.88)	(1,4-5,22,38, 42)																					
A	US,A, 4434236 (FREYTAG) 28 February 1984 (28.02.84)	(1,5)																					
A	WO,A, 84/02193 (E.I. DUPONT DE NEMOURS AND COMPANY) (EBERSOLE et al) 7 June 1984 (07.06.84)	(1-2,4-5)																					
A	AU,A, 66596/86 (ORTHODIAGNOSTIC SYSTEMS INC)(FINKENAU) 18 June 1987 (18.06.87)	(1-2,4-5,13, 38)																					
A	EP,A, 0201211 (WHITTAKER CORPORATION)(KHEILL et al) 12 November 1986 (12.11.86)	(1,4-6,22)																					
* Special categories of cited documents: ** -A- document defining the general state of the art which is not considered to be of particular relevance -E- earlier document but published on or after the international filing date -L- document which may throw doubt on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified) -O- document referring to an oral disclosure, use, exhibition or other means -P- document published prior to the international filing date but later than the priority date claimed -T- later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention -X- document of particular relevance: the claimed invention cannot be considered novel or cannot be considered to involve an inventive step -Y- document of particular relevance: the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art. -A- document member of the same patent family																							
<b>IV. CERTIFICATION</b> <table border="1"> <tr> <td>Date of the Actual Completion of the International Search 10 November 1988 (10.11.88)</td> <td>Date of Making of this International Search Report 17 NOVEMBER 1988 (17.11.88)</td> </tr> <tr> <td>International Searching Authority Australian Patent Office</td> <td>Signature of Authorized Officer  J.H. CHAN</td> </tr> </table>			Date of the Actual Completion of the International Search 10 November 1988 (10.11.88)	Date of Making of this International Search Report 17 NOVEMBER 1988 (17.11.88)	International Searching Authority Australian Patent Office	Signature of Authorized Officer  J.H. CHAN																	
Date of the Actual Completion of the International Search 10 November 1988 (10.11.88)	Date of Making of this International Search Report 17 NOVEMBER 1988 (17.11.88)																						
International Searching Authority Australian Patent Office	Signature of Authorized Officer  J.H. CHAN																						

ANNEX TO THE INTERNATIONAL SEARCH REPORT ON  
INTERNATIONAL APPLICATION NO. PCT/AU 88/00274

This Annex lists the known "A" publication level patent family members relating to the patent documents cited in the above-mentioned international search report. The Australian Patent Office is in no way liable for these particulars which are merely given for the purpose of information.

Patent Document Cited in Search Report		Patent Family Members			
US	4752562	CA 1179940	EP 70300	ES 508875	
		ES 8301367	IT 1147804	US 4517288	
		US 4774174	WO 8202601		
AU	53187/86	WO 8604352	EP 209554		
WO	8402193	DK 3759/84	EP 126772		
EP	201211	JP 61286755			
AU	66596/86	EP 228225	JP 62188971		

END OF ANNEX



**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

☐ BLACK BORDERS

☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES

☐ FADED TEXT OR DRAWING

☒ BLURRED OR ILLEGIBLE TEXT OR DRAWING

☐ SKEWED/SLANTED IMAGES

☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS

☐ GRAY SCALE DOCUMENTS

☐ LINES OR MARKS ON ORIGINAL DOCUMENT

☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY

☐ OTHER: \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**